



READ THE MANUAL BEFORE OPERATING THIS PRODUCT


COMPEX  [®]

MOLECULE[™]

Vibrating Massage Ball


INTRODUCTION

Pinpoint Pressure - The Complex Molecule vibrating massage ball has four levels of vibration intensity to relieve pain from stiff or aching muscles. Increases circulation in targeted areas, relieve soreness, and inflation to muscles. Use this recovery tool before and after a workout to heal and nourish your body.

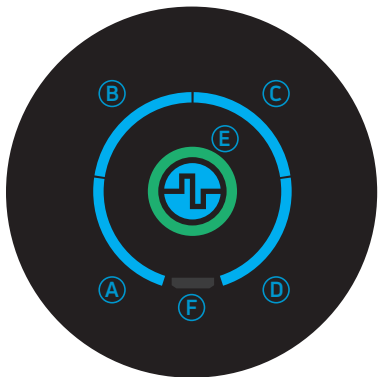

- (A)** Level 1 - 1700rpm
- (B)** Level 2 - 1900rpm
- (C)** Level 3 - 2200rpm
- (D)** Level 4 - 2200-1700rpm
- (E)** Power Indicator
- (F)** USB charging port (5V 2A)

POWER INDICATOR

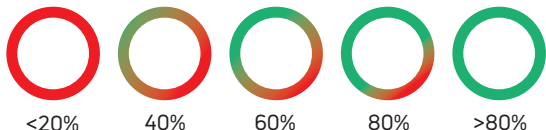

BASIC SPECIFICATION		SPEED 1	SPEED 2	SPEED 3	SPEED 4	REMARK
Noise(dB)		<54db	-	-	54db	Condition: 30cm distance
Rotational speed (with loading)		1700	1900	2200	1700-2200	
Hertz (Hz)		28.33 Hz	31.67 Hz	36.67 Hz	28.3Hz-36.6Hz	
Power consumption testing	Voltage (V)	8.1	8.06	7.88	7.88	
	Current (A)	0.29	0.41	0.91	0.8	
	Power (W)	2.35W	3.3W	7.2W	6.3W	
Charge and Discharge test	Charge time	1.5H				
	Discharge time				1-1.5H	Every 15 minutes of work automatically stops for one round
Stand-by Power Consumption		<60uA				
Adapter	AC100~240V 50/60Hz, 5V 2A US standard USB cable (Certification: UL/ETL ROHS)					
Battery	7.4V 1200mAh Li-Polymer Battery (2 in series)					

FUNCTION INSTRUCTIONS

The product is designed to have four levels of vibration intensity. To power on/off, press down the power button for 3 seconds until the power indicator light illuminate/turn-off. To select the desirable vibration level, press the power button to cycle from one level to the next. The blue light will indicate the vibration level you are on.

CAUTION

- It is not recommended to deplete the battery less than 20% for a long period of time. Please fully charge the device to avoid discharge of the battery when the device is lower than 25% battery life.
- In each mode, the device will automatically turn off after 15 minutes of use.

SAFETY WARNING

THIS PRODUCT MUST BE USED IN ACCORDANCE WITH THE FOLLOWING INSTRUCTIONS:

- Used by and/or administered by adults only (21 and over).
- Stop immediately if you feel any pain or discomfort.
- Avoid contact with head and/or other bony areas.
- Use on body parts for no more than 1 minute per specific location. Keep moving over muscles and avoid bony areas.
- Avoid putting fingers, hands and/or hair near the ball/shaft section of Molecule™.
- Do not drop Molecule™.
- This appliance is not intended for use by persons (including

children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.

- Children should be supervised to ensure that they do not play with the appliance.
- The appliance is only to be used with the power supply unit<model ICP15A-050-2000B> provided with the appliance.
- WARNING: For the purposes of recharging the battery, only use the detachable supply unit provided with this appliance.

CONTRAINDICATIONS

PLEASE DO NOT USE MOLECULE™ WITHOUT RECEIVING APPROVAL FROM YOUR PHYSICIAN IF THE FOLLOWING CONDITIONS APPLY:

- Recent surgeries of any kind
- Skin conditions and/or wounds affected by vibrating massagers
- Herniated disks
- Pregnancy
- Diabetes/neuropathy or retinal damage
- Wear of pacemakers
- Injuries affected by vibrating massagers
- Epilepsy or migraines
- Disc herniation
- Spondylolisthesis
- Vertebral detachment or spine joint rigidity
- Recent replacement of joints or IUD
- Metal pins or metal plates
- Any health-related diseases

WARNINGS

- Do not use if injured.
- Do not immerse in water, sand and/or loose debris that might get into motor housing vent.
- Do not remove screws or attempt to disassemble.
- This appliance contains batteries that are non-replaceable

ADVERSE REACTIONS

You may experience skin irritation and mild bruising if Molecule™ is used in one spot for more than 60 seconds or over bony areas. Adverse reactions may be experienced by people recovering from surgery, skin conditions, wounds, nerve conditions, broken bones and/or contusions. To avoid reactions, users should consult doctor prior to using Molecule™.

CHARGING INSTRUCTIONS

1. Please fully charge the device before first use.
2. While charging, the power indicator will turn to green according to the actual power level of the battery.
3. The product cannot be used while charging.
4. Please only use the adapter provided. Do not use other adapters to avoid damage to the battery and potential danger.
5. Charging time: 1.5 hours.

DECLARATION OF CONFORMITY

FCC PART 15B Compliance Statement

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

COMPEX VIBRATING MASSAGE BALL

MODEL: CX202WL03

Powered by Lithium-polymer battery (7.4V , 1200mAh, 8.88Wh)

Rated Input: 5V 2A

Rated Power: Max. 12W

Use only with <ICP15A-050-2000B> supply unit



Made in China

WARRANTY

This device carries a limited warranty of one year from the date of delivery. The warranty applies to the device only, and the accessories are not covered by this warranty.

During the warranty period, defective items will be repaired or replaced at no charge. Any evidence of misuse, abuse, alterations, or externally caused damage may have this warranty invalid.

For more information, please contact the distributor/manufacturer.

CONTACT INFORMATION

Distributed By DJO, LLC
1430 Decision Street • Vista, CA 92081
1-877-266-7398 (877-COMPEX8)

Copyright © 2020 by DJO, LLC
Individual results may vary. Neither DJO, LLC nor any of its subsidiaries dispense medical advice. The contents of this document do not constitute medical advice. Rather, please consult your healthcare professional for information on the courses of treatment, if any, which may be appropriate for you.

DJO, LLC
A DJO Global Company
1430 Decision Street
Vista, CA 92081-8553
U.S.A. | [complex.com](https://www.complex.com)

CX2021F09 Rev A